INSTRUCTIONS FOR USE – EdgeOne GlidePath™

EN	INSTRUCTIONS FOR USE	2
BG	ИНСТРУКЦИИ ЗА УПОТРЕБА.....	5
DA	BRUGSANVISNING	8
DE	GEBRAUCHSANWEISUNG	11
EL	ΟΔΗΓΙΕΣ ΧΡΗΣΗΣ	14
ES	INSTRUCCIONES DE USO.....	17
ET	KASUTUSJUHEND	20
FI	KÄYTTÖOHJE	23
FR	NOTICE D’UTILISATION	26
HR	UPUTE ZA UPORABU	29
HU	HASZNÁLATI UTASÍTÁS	32
IT	ISTRUZIONI PER L’USO.....	35
LT	NAUDOJIMO INSTRUKCIJA.....	38
LV	LIETOŠANAS INSTRUKCIJAS.....	41
NL	GEBRUIKSAANWIJZING	44
NO	BRUKSANVISNING	47
PL	INSTRUKCJA UŻYWANIA	50
PT	INSTRUÇÕES DE UTILIZAÇÃO	53
RO	INSTRUCȚIUNI DE UTILIZARE.....	56
SK	NÁVOD NA POUŽITIE	59
SL	NAVODILA ZA UPORABO	62
SV	BRUKSANVISNING	65
TR	KULLANIM TALİMATLARI	68


FOR DENTAL USE ONLY

INSTRUCTIONS FOR USE – EdgeOne GlidePath™

1) GENERAL DESCRIPTION

Composition:

The instrument is made of nickel-titanium alloy and the endo stop is made of silicone.

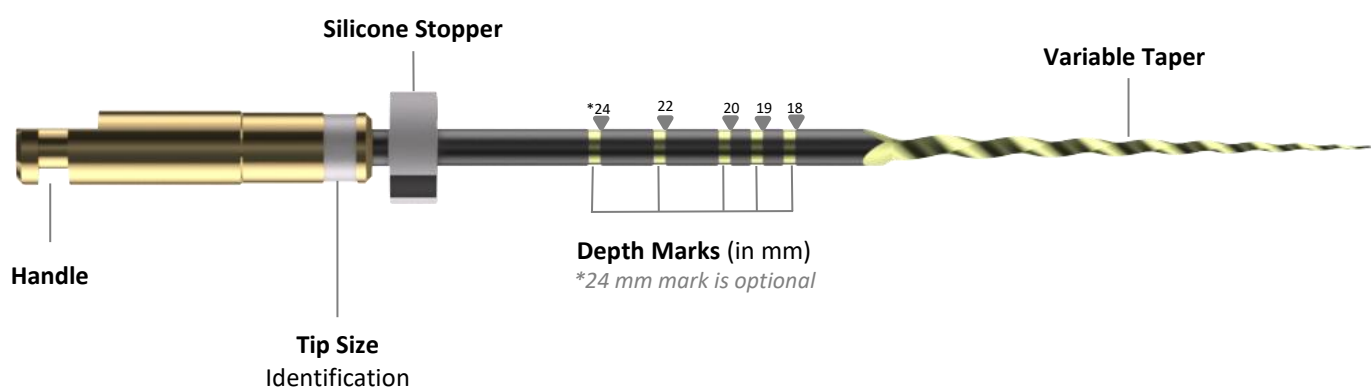The tip size is indicated by the color band. The taper is variable.

These instruments are compatible with endodontic motors compliant with ISO 1797.

Recommended speed: 350 rpm in counterclockwise motion; torque limit: 300 g·cm (2.94 N·cm).

Size	Taper	Length
016	.02v	21mm
016	.02v	25mm
016	.02v	31mm

2) INTENDED USE

The EdgeOne GlidePath™ endodontic instruments are intended for root canal treatment.

3) INDICATIONS FOR USE

The EdgeOne GlidePath™ are indicated in case of pulp or periapical diseases that require root canal treatment. All instruments are designed for preparation and shaping of root canals.

4) CONTRAINDICATIONS AND ADVERSE REACTIONS

None known.

5) WARNINGS AND PRECAUTIONS

- These instruments have not been tested on children, pregnant nor breastfeeding women.
- These endodontic instruments must not be used on individuals with a known allergic sensitivity to nickel-titanium and silicone.
- Single use instruments must not be reprocessed for subsequent use nor be reused. This could compromise their integrity.
- Each instrument may be used to treat up to four (4) canals during a single procedure. Discard after use.
- For your own safety, wear personal protective equipment (masks, gloves, and goggles) during the treatment and reprocessing of the instrument.
- Do not use if sterile packaging seal is broken, damaged or wet.
- Do not use the instruments after expiration date.
- Reprocess the instrument appropriately before disposal.
- Carefully read the label, marking on the packaging and this instruction for use to ensure the correct identification and use of device.
- Use a dental dam when using endodontic instruments to avoid inhalation or ingestion by the patient.
- Take several radiographs from different angles to determine the anatomy of the root canals (length, width and curvature).
- Always inspect the instrument(s) before use, and discard (it) them if there are any visible defect(s).
- Cracks, deformations, signs of corrosion, loss in color or marking are signs that the device is no longer able to achieve the required performance level and should be discarded.
- Regularly clean flutes and inspect instrument during use, and discard if it shows any signs of distortion or wear.
- The packaging of the instruments should be carefully examined before opening (packaging integrity, no humidity, no visual degradation of instruments) to ensure that the packaging's integrity has not been compromised during transport and storage.
- If the packaging has been opened, damaged, become wet or the expiry date has passed, the sterile state of the instruments inside the packaging is not guaranteed. Discard the device(s) that are single use and consult the reprocessing instructions available on EdgeEndo's website for multiple use device(s).
- All serious events occurring in relation to the product must be reported to the manufacturer and the competent authority according to local regulations.

6) PERFORMANCE CHARACTERISTICS

The EdgeOne GlidePath™ files are compliant with all applicable testing requirements of the ISO 3630 series of standards.

INSTRUCTIONS FOR USE

EN

7) PROTOCOL FOR USE

1) Straight Line Access:

- Form a straightline access to the canal orifice.
- Use a #10 K-file to establish the working length with an apex locator and/or radiograph.
- Confirm patency with the #10 K-File.
- Irrigate frequently along the process.
- Use The EdgeOne GlidePath™ in one or more passes until the working length is reached.
- Shape the canal with EdgeEndo rotary or reciprocating files.

2) Final Irrigation and Drying:

Perform a final irrigation protocol according to clinical preference. Dry the canal with paper points that correspond to the final file size taken to working length.

3) Obturation:

Obturate the canal system using a matching gutta-percha cone or thermal carrier system.

8) CONDITIONS AND TERM OF STORAGE, SHELF-LIFE

The shelf-life is 5 years for sterile instruments. The use-by date is indicated on the packaging. Sterility is preserved and guaranteed for the whole shelf life unless the packaging is broken, or the storage conditions are not respected.

9) DISPOSAL

When a device reaches the end of its life, make sure that it is discarded in accordance with the applicable laws and regulations.

10) SYMBOLS

For explanation of symbols for IFUs and labels, please consult the document “Glossary of the symbols used by EdgeEndo” available on EdgeEndo’s website.

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Henry Schein Services GmbH,